


ALUMINIUM OVERBED FOLDABLE TABLE


UNFOLD




01 Place case on a flat surface and open cap.

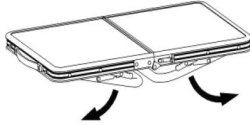

02 Unfold case & pull legs out.

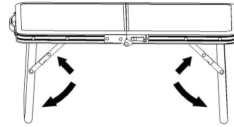


03 Straighten legs & lock using safety lock.

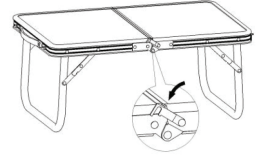


04 Close lock on both sides to prevent folding.

FOLDING

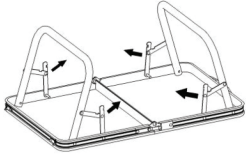

01 Open the safety lock on legs of table

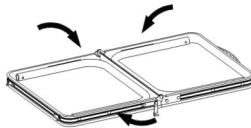

02 Fold both legs inwards



03 Unlock centre lock & Close table.



04 Close cap. You are ready to go.

CAUTION



Please use on flat surface .



Please do not sit on table.



Take care of hand when you open or close.



Please do not expose to extreme moisture or heat.

Import & Distributed by



ETERNITY
DISTRIBUTION PVT LTD

No. 154, 1/F, 3rd Main Road,
Palaniappa Nagar, Pudur, Ambattur,
Chennai - 600 053. Tamil Nadu, INDIA.

+91 44 2686 3690

www.ernitydistribution.com

marketing@ernitydistribution.com